



Octet QKe System

Enhanced quantitation and kinetics performance with biosensor re-racking

Key features

- Increased sensitivity and dynamic range for protein quantitation
- Excellent signal resolution for measuring protein and peptide kinetics
- Greater operational cost savings with biosensor re-racking and regeneration

The Octet® QKe system is the enhanced-performance version of the first-generation Octet QK instrument from the Octet family of label-free biomolecular interaction analysis instruments. Two acquisition rate settings and a higher-performance spectrometer combine to provide greater sensitivity and a wider dynamic range for protein quantitation and kinetic analysis. In addition, biosensor re-racking provides enhanced assay flexibility and operational cost savings.

Increased sensitivity and dynamic range for protein quantitation and kinetics

The Octet QKe instrument contains the same high-performance spectrometer used in the Octet RED96 system. In addition, the Octet QKe system provides a High Sensitivity acquisition

rate setting for both quantitation and kinetics assays, allowing measurement of lower ranges of protein concentrations and smaller molecular weight analytes than was possible with the Octet QK system.

The Standard acquisition rate of 0.6 Hz with 5 averages per data point (1 data point per 1.6 seconds) is the default setting on the Octet QKe system and equivalent to that on the Octet QK system. The High Sensitivity acquisition rate of 0.3 Hz (1 data point per 3.3 seconds) allows the software to perform 40 averages per data point to reduce noise and enhance signal to noise ratios, shown in Table 1, and gives the Octet QKe system enhanced sensitivity for quantitation of proteins and peptides.

Kinetic analysis of proteins and peptides benefits from these Octet QKe system enhancements over the Octet QK system in two ways:

- Lower molecular weight proteins and peptides can be detected and analyzed on the Octet QKe system.
- The binding data for proteins or peptides at different concentrations can be better resolved on the Octet QKe system, providing reliable measurement of kinetic constants.

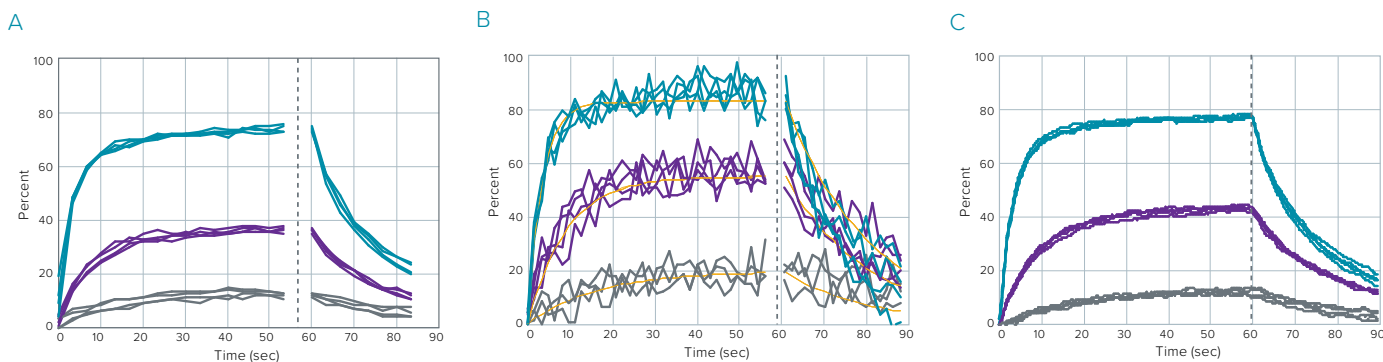


Figure 1: Comparison of the binding of insulin (5.8 kDa) to biotinylated anti-insulin antibody anchored on Super Streptavidin biosensors on the Octet QKe system (A), the Octet QK system (B), and the Octet RED system (C). The Octet QKe system data was obtained using the Low acquisition rate setting of 0.3 Hz and 40 averages. The Octet RED system data was obtained using 5 Hz acquisition rate and 20 averages and the Octet QK system data was obtained using a 0.6 Hz acquisition rate and 5 averages.

System	Octet QK	Octet QK ^e	Octet RED
Setting	0.6 Hz, 5 avg	0.3 Hz, 40 avg	0.6 Hz, 5 avg 5 Hz, 20 avg
k_{on} (1/Ms)	1.51E+05	1.37E+05	1.08E+05 1.03E+05
k_{off} (1/s)	7.48E-02	8.00E-02	7.25E-02 6.56E-02
K_D (M)	4.95E-07	5.82E-07	6.71E-07 6.39E-07
S/N @ [0.5 K_D]	4.7	30.8	14.3 29.1

Table 1: Comparison of performance of the Octet QK, QK^e (at two different acquisition rates) and RED systems for the insulin (5.8 kDa)-anti insulin interaction.

Octet QK^e system specifications*

Sample and analysis	
Detection Technology	Bio-Layer Interferometry (BLI)
Biosensor Type	Disposable, single-use fiber optic biosensors with optional reuse by regeneration and/or re-racking
Information Provided	<ul style="list-style-type: none"> Kinetic and affinity analysis (k_{obs}, k_a, k_d, K_D) Concentration data analysis (no need for background subtraction) Epitope binning and cross-blocking matrices and trace overlays
Data Presentation	<ul style="list-style-type: none"> Plots displaying kinetic binding, equation fits, and residuals of fits Tabulated kinetic and quantitation data
Sample Types	Proteins, antibodies, peptides, media containing serum, buffers containing DMSO, periplasmic fractions, untreated cell culture supernatants, and crude cell lysates
Sample Plate	Standard 96-well, black, flat bottom microplate
Sample Volume	<ul style="list-style-type: none"> 180–220 μL/well (96-well microplate) Nondestructive testing, easily recovered
Orbital Flow Capacity	Static or 100–1500 rpm
Analysis Temperature Range	(Ambient + 4°C)–40°C, 1°C increments

*All specifications are subject to change without notice.

Distributore per l'Italia:

ALFATEST
strumentazione scientifica

alfatest@alfatest.it
06 87465556/7

www.alfatest.it
www.alfatestbio.it
www.alfatestlab.com

FORTÉBIO

www.fortebio.com

ForteBio
47661 Fremont Boulevard
Fremont, CA 94538
888.OCTET-75 or 650.322.1360
fortebio.info@moldev.com

ForteBio Analytics (Shanghai) Co., Ltd.
No. 88 Shang Ke Road
Zhangjiang Hi-tech Park
Shanghai, China 201210
salesops.china@moldev.com

Molecular Devices (UK) Ltd.
660-665 Eskdale
Winnersh Triangle
Wokingham, Berkshire
RG41 5TS, United Kingdom
+44 118 944 8000
uk@moldev.com

Molecular Devices (Germany) GmbH
Bismarckring 39
88400 Biberach an der Riss
Germany
+00800 665 32860

In comparison, the Octet QK384 system provides equivalent performance at greater throughput. The Octet HTX, RED96 and RED384 systems provide the best performance among the Octet family of instruments for quantitation and kinetic analysis, with the broadest dynamic range, greatest sensitivity, and ability to measure small molecule: protein binding and fast binding interactions.

Making quality analysis affordable

The Octet QK^e system provides re-racking and reuse of biosensors at the end of an experiment. Biosensor re-racking provides enhanced flexibility for loading ligands on biosensors and reduces operational costs by regenerating and re-using biosensors. Combined with superior data quality, the Octet QK^e system is a value-added alternative to other label-free instruments for your assays.

Quantitation and kinetics

Throughput Up to 8 assays in parallel; up to 96 assays per 96-well microplate

Analysis time per sample

- hIgG quantitation in 2 minutes for 8 samples, \leq 32 minutes for 96 samples
- Real-time kinetic analysis experiments from 5 minutes to 4 hours

Baseline noise

- \leq 3 pm (RMS) at High Sensitivity acquisition rate
- \leq 8 pm (RMS) at Standard acquisition rate

Quantitation range for hIgG

- 1 μ g/mL to 700 μ g/mL at low rpm
LOD of 0.05 μ g/mL at 1000 rpm

Physical specs

Dimensions 18.6 in (H) x 17 in (D) x 20.8 in (W)
47 cm (H) x 43 cm (D) x 53 cm (W)

Weight 54 lb (24.5 kg)

Electrical requirements

- Mains: AC 100–240 V, 5.0–2.0 A, 50/60 Hz, single phase
- Power consumption: 120 W (240 W peak)

Safety standards

- CE, Nemko

Ordering information

Part No.	UOM	Description
OCTET QKE	System	Includes Octet QK ^e instrument, Octet software, desktop computer, LCD monitor, accessory kit, and one-year warranty